



## Conditions

Column: ACE 3 C18  
Dimensions: 50 x 4.6 mm  
Part Number: ACE-111-0546  
Mobile Phase: H<sub>2</sub>O/MeCN/TFA (64:34:0.5 v/v/v)  
Flow Rate: 2 mL/min  
Injection: 20 µL  
Temperature: 28 °C  
Detection: UV, 254 nm

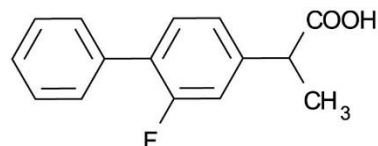

Flurbiprofen

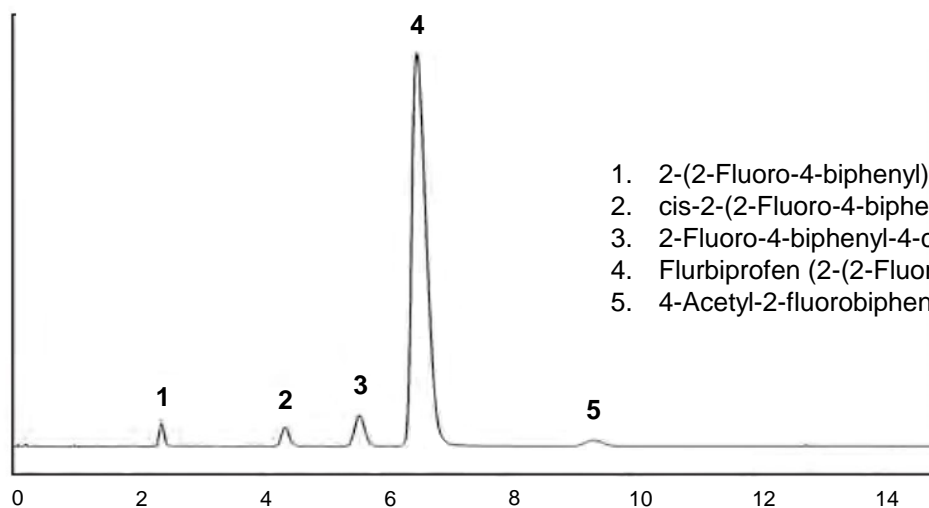

1. 2-(2-Fluoro-4-biphenyl)-2-hydroxypropionic acid
2. cis-2-(2-Fluoro-4-biphenyl)-2-hydroxypropionic acid
3. 2-Fluoro-4-biphenyl-4-carboxylic acid
4. Flurbiprofen (2-(2-Fluoro-4-biphenyl)propionic acid)
5. 4-Acetyl-2-fluorobiphenyl