

Rapid analysis of Cold and Flu Medicine Components

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN4780

Conditions

Column: ACE Excel 2 C18-AR
Dimensions: 100 x 3.0 mm
Part Number: EXL-109-1003U
Mobile Phase: A: 0.1% formic acid in H₂O
B: 0.1% formic acid in MeCN

Time (mins)	%B
0	5
3.6	20
5.4	80
6	80
6.5	5

Flow Rate: 1.0 mL/min
Injection: 5 µL
Temperature: 40 °C
Detection: UV, 275 nm
Instrument: Chromaster UltraRs

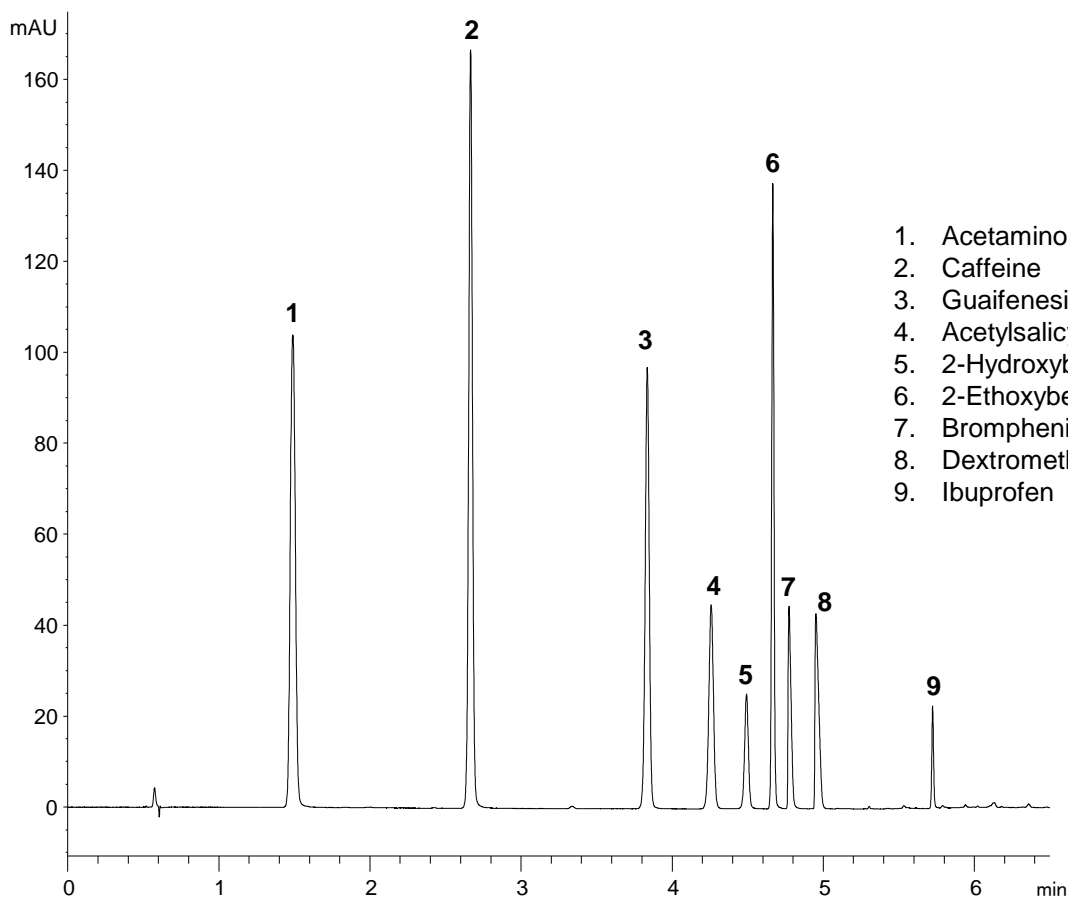

© 2017 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com

