


USP Method for Analysis of Povidone Impurity B (2-Pyrrolidone)

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN6840

Conditions


Columns:	ACE 5 C18
Dimensions:	150 x 4.6 mm
Part Number:	ACE-121-1546
Mobile Phase:	5:95 MeOH/H ₂ O
Flow Rate:	0.8 mL/min
Injection:	5 µL
Temperature:	40 °C
Detection:	UV, 205 nm

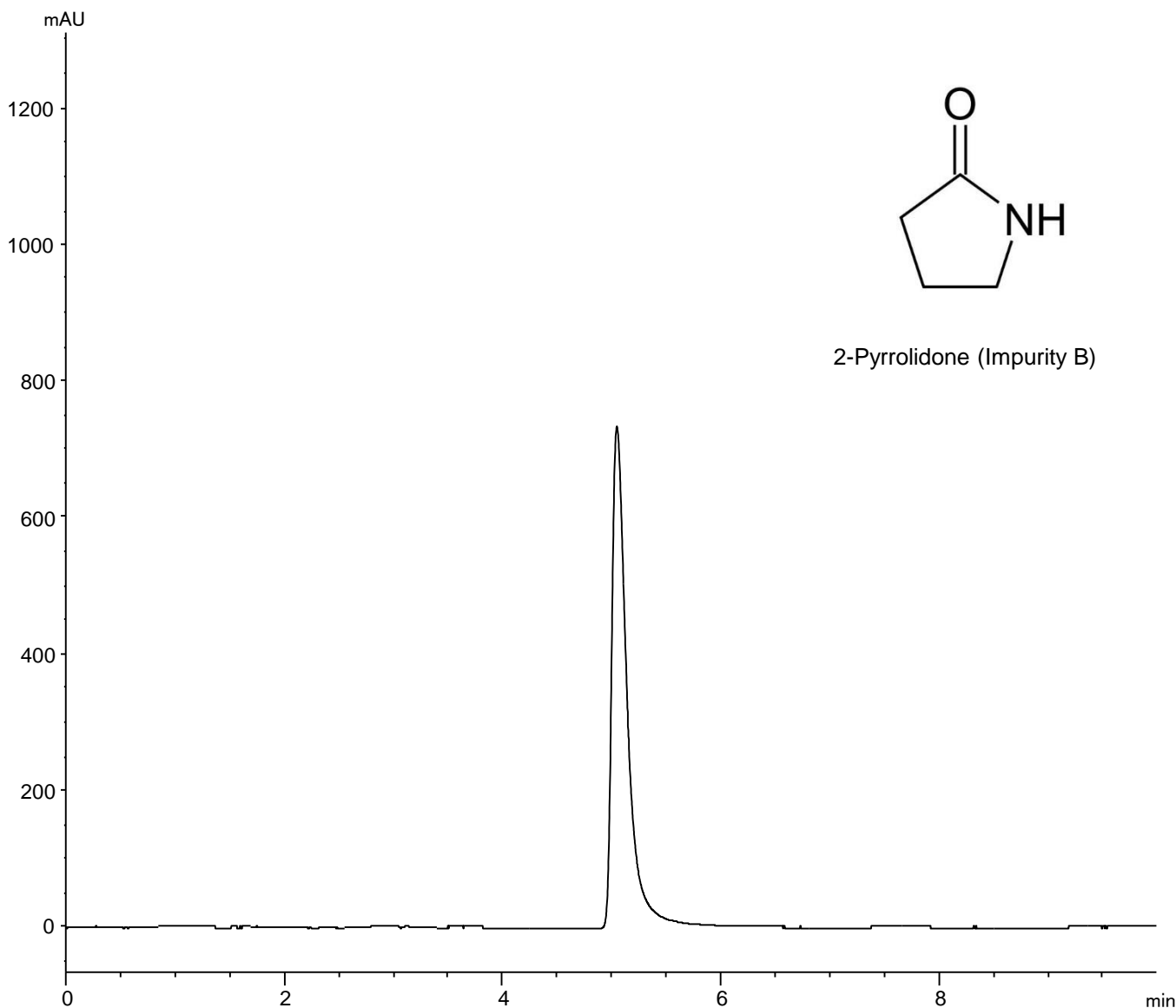

Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com



ACE[®]
UHPLC & HPLC
Columns